




## XLR Plus

The XLR Plus is a mid-size power chair that is design for the daily user. It sports a mid-wheel drive system and has the added stability of having six points of contact with the ground. The XLR Plus also has our stable four post seating system and full independent suspension system for smooth drive operation that will handle most driving surfaces. The unit has easy battery access, adjustable footplate, arm rest and 300 lb. weight capacity.

- + 300lb Weight Capacity
- + Independent Front & Rear Suspension
- + Four Post Seating System
- + Exceeds Medicare's Specs for Indoor/Outdoor use
- + Reclining Captain Seat

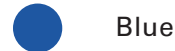

Model Number: 858WM  
Medicare Billing Code: K0823

HCPCS codes provided should not be considered as legal advice, and does not guarantee reimbursement.

## Technical Specifications

Overall Dimensions (L x W x H)	in	36" x 22" x 38.6"
Number, Size of Tires	Front, Rear	2, 6" (2, 10") 2, 5"
Suggested User Weight <sup>1</sup>	lb	Max. 300
Weight of Heaviest Part <sup>2</sup>	lb	-
Weight <sup>3</sup>	lb	205
Each Battery Module Weight	lb	34
Battery Size	Type x PCS	12V50AH x 2
Charger	Amp	Off board (5A)
Output Power of Motor	HP	1.34 x 2
Maximum Speed	MPH	5
Maximum Climbing Angle/ Safe Climbing Angle <sup>4</sup>	%(°)	10(6)
Range (Per Charge) <sup>6</sup>	Mile	21.3 - 23.8
Turning Radius	in	18"
Model #		858WM
Standard Seat	W x D x BH	SV191821B

**Note:** <sup>1</sup> On level Road. <sup>2</sup> Without basket and seating. <sup>3</sup> With basket, seat & batteries. <sup>4</sup> Base on user weight of 200LB. <sup>5</sup> Based on user weight of 165LB. <sup>6</sup> After the batteries and mechanical moving parts fully break in.